



MediSim VR Setup Guide

Complete Guide to Setting Up Your Healthcare Training Platform

For any questions or inquiries not addressed
in this guide, please contact our support team
immediately at contact@medisimvr.com.



Contents

1 | Setup

VR Hygiene & Safety	03
Hardware Connection and Software Preparation	04
Launching MediSim VR Application	07

2 | Session Modes Explained

Guided Mode	09
Evaluation/Assessment Mode	10

3 | Tips and Tricks

Adjusting Headset Fit	12
Audio Settings	13

4 | Troubleshooting

Common Issues and Solutions	13
Support Contacts	15

VR Hygiene

Maintaining proper hygiene is essential when using shared VR equipment, especially in healthcare training environments.

Before Each Use:

- Wash hands thoroughly or use hand sanitizer
- Consider wearing nitrile gloves while cleaning equipment
- Wear a scrub cap when appropriate for your training scenario



Cleaning the Equipment:

- Use non-abrasive anti-bacterial wipes to clean all surfaces that contact users
- **IMPORTANT:** Do not use alcohol wipes on headset lenses or tracking cameras
- Use a dry microfiber cloth only for the lenses
- Allow equipment to completely air dry (minimum 10 minutes) before next use



Before Hardware Connection


Before proceeding with hardware connections, ensure the following:



- The computer is plugged into a power outlet and ready to switch on
- Monitor is connected and turned on



- VR headset (Meta Quest) is available at the station and ready to use
- Link Cable and HDMI cable are accessible

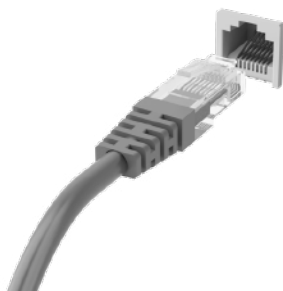

- LAN (Ethernet) cable is available for internet connection
- All USB ports are functional (preferably USB 3.0 or Type-C)
- Space is cleared for safe VR usage

Step-by-Step Hardware Setup Flow

1. Power Connection



(A)

Ensure the computer is properly connected to a power source using its power cable.



(B)

Plug the computer into a working power outlet and switch on the power (If connected to a UPS, wait until the initial beep sequence is complete before switching on the power)



(C)

The VR headset (e.g., Meta Quest) will receive power directly through the Link cable once it is connected to the PC.

Step-by-Step Hardware Setup Flow

2. Switch On the PC




(A)

Press the power button on the desktop



(B)

Wait for the system to boot up to the Windows login screen.

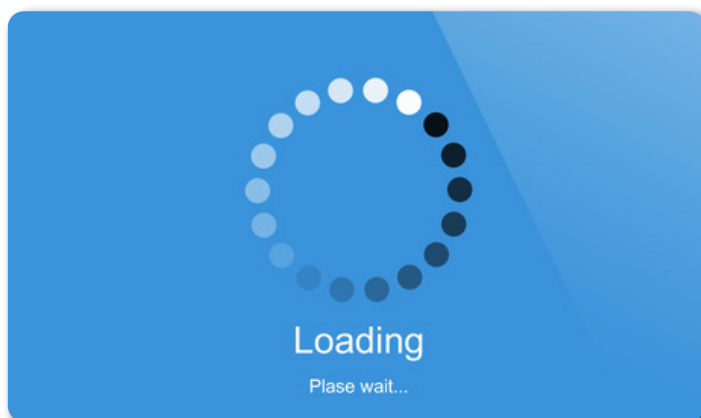

(C)

Log in with your credentials and wait for the desktop (home screen) to fully load.

Hardware Connection and Software Preparation:

1.Computer Setup

A. Power Connection:

Connect the power cables to both the VR headset and the computer.



B. Connect to Computer:

Use the Link cable to connect your VR headset to your computer.



- Plug one end into the VR headset.
- Plug the other end into your computer's USB-C or USB 3.0 port (depending on your device):

Laptop: Use the Type-C port.

Desktop: Use the USB 3.0 port.

C. Video Connection:

Ensure the HDMI cable is connected to the graphics card.



D. Network:


Connect the computer to the LAN for a reliable connection.



2. Software Setup


1. Log in to Meta Account

- If it's your first time, open a browser and go to <https://www.meta.com/setup>.
- You can create a new Meta (Oculus) account or log in to an existing one.

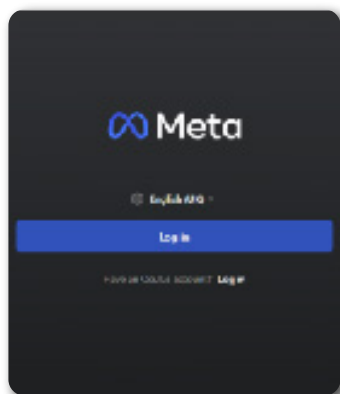


2. Install and Open the Meta Quest Link App :

- Download the **Meta Quest Link** app from the setup page.
- Once installed, **open the app** on your computer.

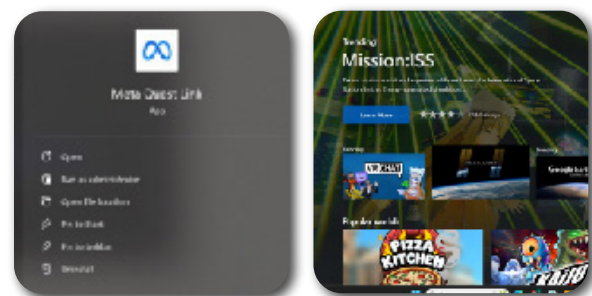


3. Linking and Configuring Your VR Headset

- Open the **Meta Quest Link** app on your computer.
- Put on your headset — a pop-up will appear asking if you want to “**Enable Link**” or “**Use Quest Link**”.
- Select “**Enable**” to switch to PC VR mode.
- You should now see the **Link home environment** (a white grid space). This confirms your headset is connected to the app.

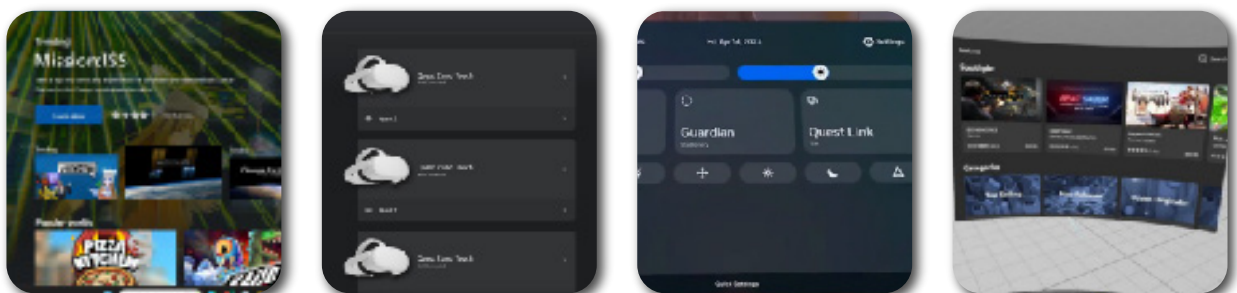



While wearing the headset, open the Settings menu, then go to Boundary.

Follow the on-screen steps to draw your play area (this ensures safe movement in VR). Make sure the boundary is clear and well-defined to avoid interruptions later.

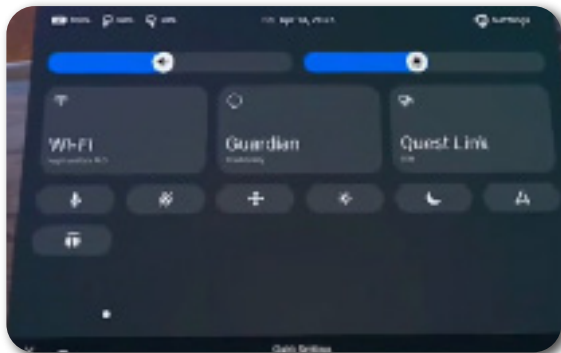

Once the boundary is set, take off the headset and look at the Meta Quest Link app on your PC.

You should now see your device listed as “Connected” under the Devices tab. This means your headset is successfully linked and ready for use.

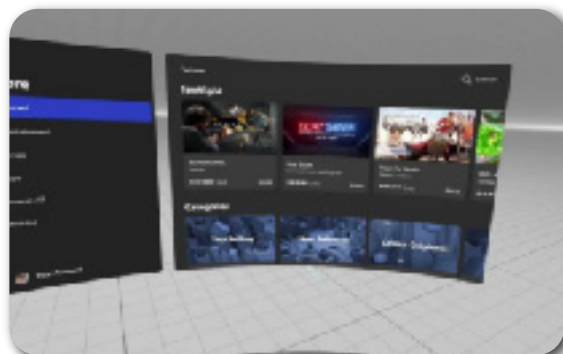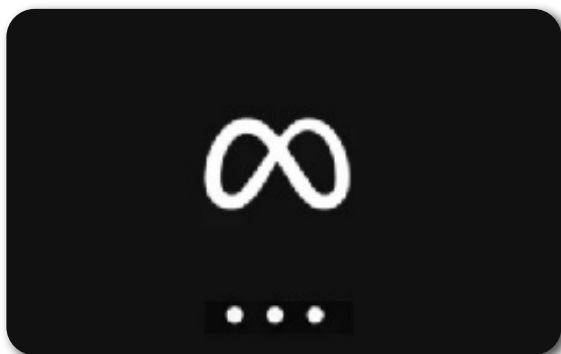

After Software Connection


After launching Quest Link, allow 3–4 minutes for the Meta app to fully load and establish a stable connection. Avoid interrupting the process during this time.

TIP

- Before starting any training session, make sure your headset is connected to Wi-Fi. Then, restart the headset to check for any system or app updates. This helps ensure a smooth, uninterrupted training experience.

Using MediSim VR Application: Launch and Training Workflow



Once your headset is connected and the Link environment is active, launch the MediSim VR application from your computer. The headset will automatically switch to the app, and you'll arrive at the login screen.

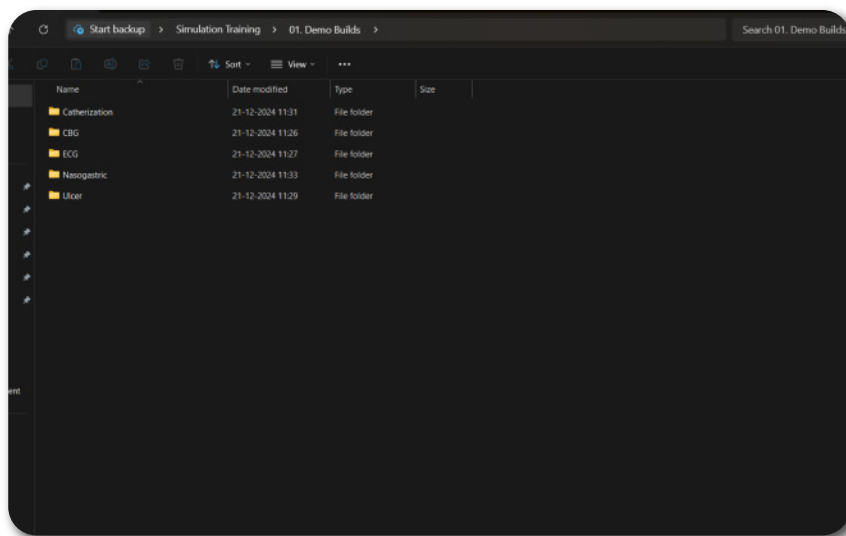

A suggested roadmap for new users and trainees.


Example Flow:

1. Start with **Guided Mode**
2. Repeat until you achieve 100% success in critical steps
3. Move to **Evaluation Mode**
4. Review performance dashboard after each session
5. Seek feedback from instructors (if applicable) t

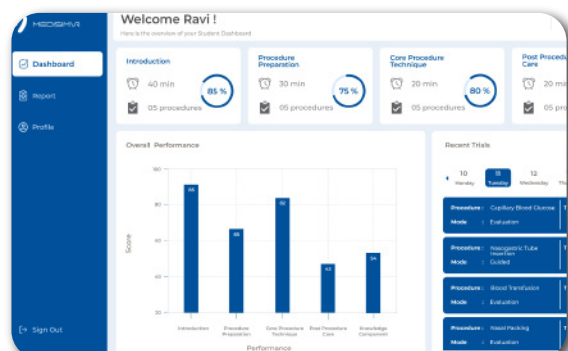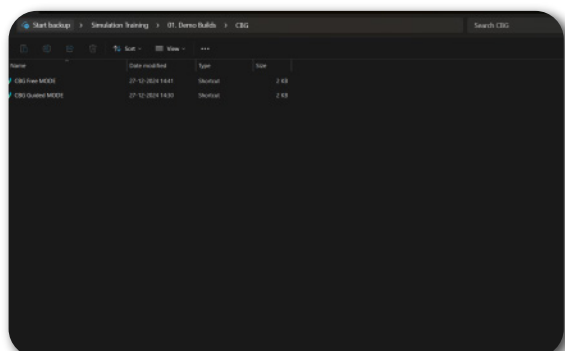


Connecting to MediSim VR

1. With the headset on and controllers in hand, click to launch the MediSim VR's Guided mode module
2. Make sure you are connected to Wi-Fi
3. Clicking the module will lead you to a log in page:

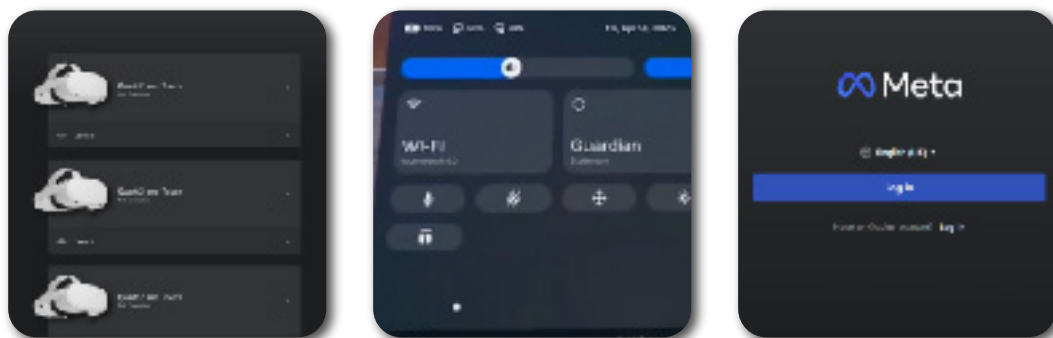

Enter your credentials and click "Log In."

MediSim VR offers two distinct learning modes designed to support users at every stage of clinical skill development:

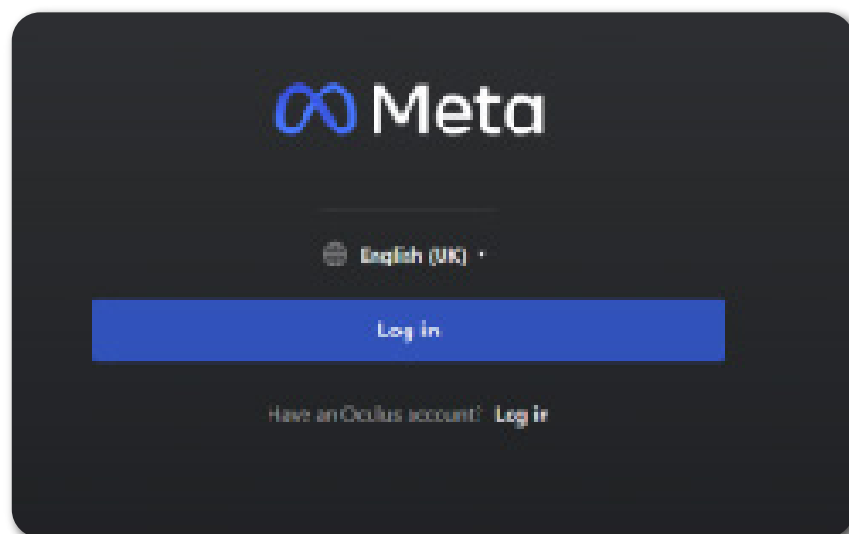

Guided Mode — Learn with Support

Guided Mode is designed to provide step-by-step instructions, visual cues, and theoretical guidance for users learning or reviewing a clinical procedure. It's ideal for first-time learners or for refreshing your knowledge.


How to Use Guided Mode

On the MediSim VR home screen, select the clinical module you want to train with.

Log in using your provided credentials.

Follow the on-screen guidance to complete each procedural step.

Once the procedure is completed, you'll be able to review your performance metrics and feedback.

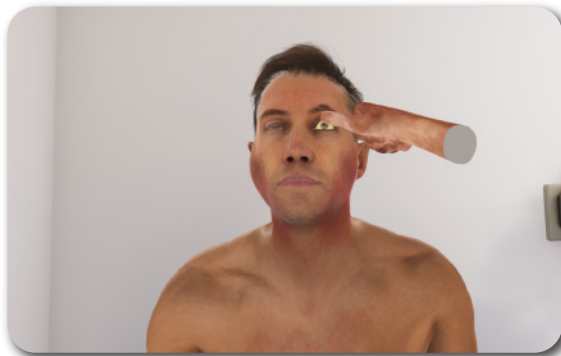

TIP

- Use Guided Mode when learning a new procedure for the first time or when you want to review and refresh your skills.

Evaluation Mode - Perform with Precision

Evaluation Mode evaluates procedural proficiency by tracking metrics like completion time and accuracy. This mode provides no guidance and is designed to test your proficiency in procedures.

How to Use Evaluation Mode

1. From the MediSim VR home screen, select the clinical module you want to be assessed on.
2. Choose the “Evaluation Mode” option.
3. Perform the procedure independently, without any step-by-step guidance.
4. Upon completion, review your detailed performance metrics, including success rate, timing, and critical step accuracy.



TIP


- Your results in Evaluation Mode are automatically saved and can be accessed anytime from your Performance Dashboard. Use this to monitor your progress and identify areas for improvement over time.

User Path

A suggested roadmap for new users and trainees.

Example Flow:

1. Start with **Guided Mode**
2. Repeat until you achieve 100% success in critical steps
3. Move to **Evaluation Mode**
4. Review performance dashboard after each session
5. Seek feedback from instructors (if applicable)

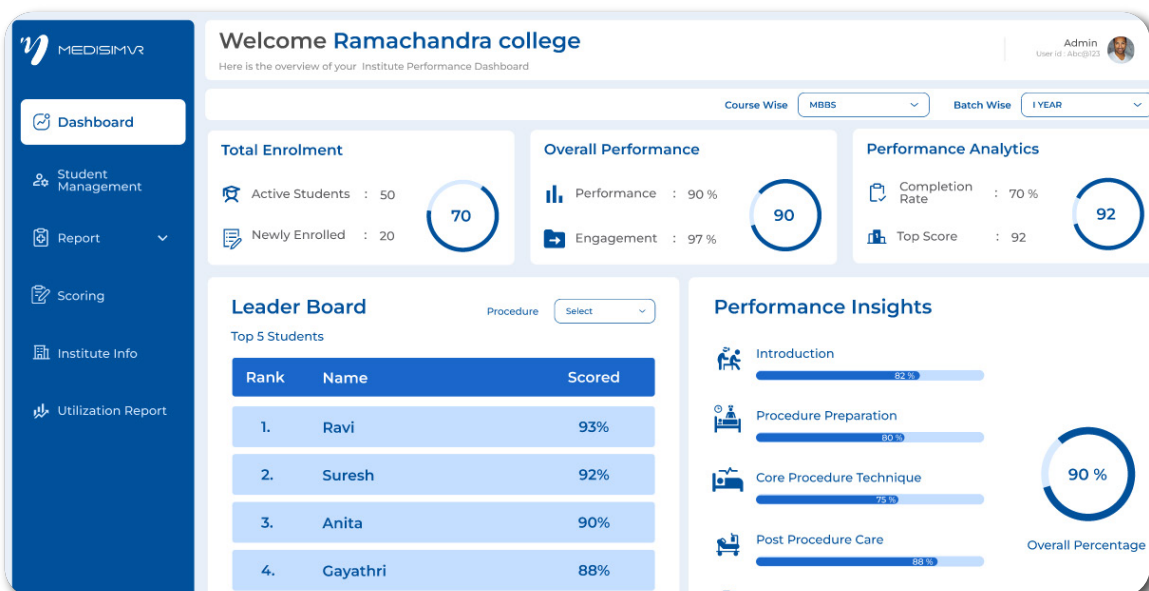

AFTER USER PATH

System Shutdown checklist

After completing your VR session, follow this checklist to properly shut down and reset the system:



Exit all VR modules by pressing Alt + F4 on the keyboard and return to the Meta Home environment



Close the Medisim VR folder with modules, opened during the session



Turn off the headset by holding the power button until the screen turns off



Ensure the headset lenses are cleaned on all surfaces that contact users



Gently return the headset to its position -ensure it's secure and not swinging freely



Close the Meta/Oculus app on the PC



Shut down the computer using the Start menu



Turn off the monitor (if applicable)



Make sure all cables are untangled and positioned safely

Practical Tips & Techniques for Comfortable VR Use and Clear Audio

For Optimal Comfort and Visual Clarity

1. Adjust the lenses until you find the best position for screen clarity

- Shift lenses right and left until they click into position
- Note the position number between the lenses for future reference

2. Adjust the side straps for a comfortable fit

- To tighten: move sliders away from the top strap
- To loosen: move sliders toward the top strap

3. Center the top strap between the two slides and on top of your head

- Adjust the Velcro for a gently snug fit

4. For users who wear glasses:

- Put the headset on your face first
- Pull the back of the headset down and tighten until comfortable




TIP

The firm arms on the headset can be moved up or down to adjust the angle against your face.

Audio Settings

Adjusting Volume

- Control the volume using the buttons on the bottom right side of the headset
- When multiple headsets are being used in the same room, turn down the volume to avoid echo effects

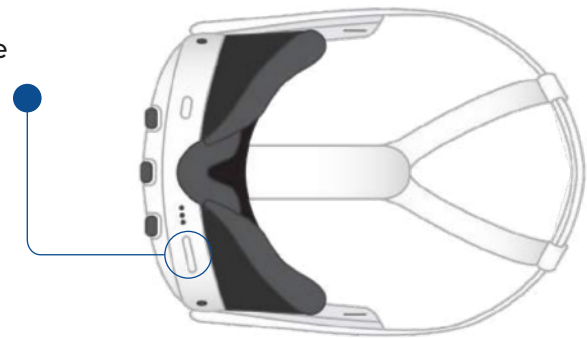

Microphone Controls

- Access the Settings menu and select 'Microphone' to mute or unmute.

Common Issues and Solutions

CHALLENGE	TROUBLESHOOTING AND SOLUTIONS
Audio Issues	Restart the headset Check volume settings on headset and PC Inspect audio output settings on PC
Login Issues	Check internet connection and Quest cable Verify internet connection Reconfirm username and password Make note of any error messages for technical support
Station Unit Not Powering On	Check power connections and cables Check power outlet functionality. Try a different power outlet or cable if possible.
Hand Controller Connectivity Issues	Check controller connectivity in device settings Verify sync between headset and controllers Check battery level (replace if below 30%) Restart the headset and reconnect the controllers if necessary

CHALLENGE	TROUBLESHOOTING AND SOLUTIONS
Link Cable Disconnects or Is Unstable	Check cable is fully plugged in both ends. Avoid sharp bends or heavy cable strain. Test with a different cable if possible. Update Oculus Link software.
PC Disconnect Errors	Reinitiate following proper instructions from the UI Update USB drivers and Oculus
PC Hardware Button Malfunction	Escalate to MediSim VR Support team
Glitches	Recheck cable connections Restart system if necessary
Graphics Issues	Escalate to support team
Software Update Failures	Ensure stable internet connection. Restart device and retry update. Free up storage space if update won't install.
Module Freezing	Check hand controls and cables Restart Try running the module on another station
User Input Errors	For accidental button presses: restart the procedure For improper tool placement: reinitiate the task following instructions
Incorrect height in VR	Reset the Guardian boundary (refer to Step 2C)
Controller Batteries Drain Quickly	Using non-rechargeable batteries Controllers staying active when not in use
Misaligned View	Hold down the Oculus Button until the circle is complete Press any button to reset the view

IMPORTANT: If challenges persist or you encounter issues not listed here, please contact our support team.

6. Safety Guidelines in VR

- Always clear the physical space of obstacles
- Use cable management to avoid tripping
- Stop immediately if you feel dizzy, nauseous, or disoriented
- Take regular breaks during long sessions

Support Contacts

For technical assistance and other inquiries, please contact:

Saurin Modh – Head of Sales
Contact: +91 9737769791
Email: saurin@medisimvr.com

General Inquiries
Email: contact@medisimvr.com
Website: [MediSim VR Website](#)
LinkedIn: [MediSim VR LinkedIn](#)

Naveen – Industry Engagement Associate
Contact: +91 9080500887
Email: naveenraja@medisimvr.com

